TOPCOM®

EAR & FOREHEAD THERMOMETER

TH-4654 + TH-4655

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This product is in compliance with the essential requirements and other relevant provisions of the R&TTE directive 1999/5/EC. This is confirmed by the CE-mark.

The Declaration of conformity can be found on: www.tristar.eu

1 INTRODUCTION

The Ear & Forehead thermometer uses infrared technology to measure infrared energy emitted from the forehead or eardrum and surrounding tissue, and then converts it into temperature value.

The Ear and Forehead Thermometer can also be used for measuring the room temperature and the temperature of food and fluids.

2 ABOUT BODY TEMPERATURE

Normal human body temperature fluctuates throughout the day; generally, lowest in the morning, highest in the evening.

Different parts of human body also have different temperatures. Therefore, the standards of measuring oral, axillary, rectal and ear temperatures are different. The human body temperature is controlled and adjusted by the hypothalamus in the brain. As the veins flow through the hypothalamus, they also flow through the ear-drum in the mid ear. Thus, measuring the ear temperature (the temperature of the ear drum membrane) can best reflect the body temperature.

Normal body temperature is a range, not a fixed value.

When using an IR ear thermometer to take temperature, the normal body temperature range is 35.8 °C - 37.6 °C. Also, people's normal temperature ranges tend to decrease with their age.

IMPORTANT NOTE

- It is important to determine your and your family members' normal temperature ranges and write them down.
- Practice to take temperatures on yourself and your healthy family members to determine their normal temperature ranges.
- It is important to clean your thermometer before the first time of use and after each application, otherwise the result may be incorrect.
- After a shower or exercise, the body temperature will rise. Take a rest before taking temperature.

3 SAFETY NOTES

WARNING


- Read carefully through the following information concerning safety and proper use. Make yourself familiar with all the functions of the equipment. Be careful to keep these advice notes and if necessary pass them on to a third party.
- This thermometer is intended for household use only. The use of this thermometer is not intended as a substitute for physician's consultation.
- The tip of the probe is the most delicate part of the thermometer, accurate readings from the thermometer can be obtained only if the tip of the probe is clean and intact.
- Don't touch or blow the infrared sensor.
- Measurement must be suspended until the ear is warmed up after outdoor activities to avoid lower results.
- After a shower or exercise, the body temperature will rise. Take a rest before taking temperature.
- For proper measurement, condition the thermometer for 30 minutes in measurement site till room temperature.
- Earwax or obstructions will affect the accuracy of the thermometer.
- It is not recommended to use the product on users suffering ear disease including otitis external and lymphantitis. Affected part is likely to deteriorate.
- Don't use the product in a wet ear hole after swimming or bathing. It is possible to damage the ear hole

DANGER


- The thermometer contains small parts (battery, battery cover) which can be swallowed by small children. For this reason, do not leave the thermometer unattended in hands of children.
- When the temperature is higher than 38°, means that the patient may have fever. Contact your doctor. Persistent fever, in particular children, has to be treated by a doctor. Get in touch with a doctor immediately.

- Keep the unit out of children's reach. A self-measurement forced by children may cause ear damage. For accidental swallowing of a battery, please contact a doctor at once.
- Stop using the product in occurrence of any pain. It is likely to damage the ear hole.
- Don't throw batteries into fire. Batteries may explode in fire

4 BUTTONS

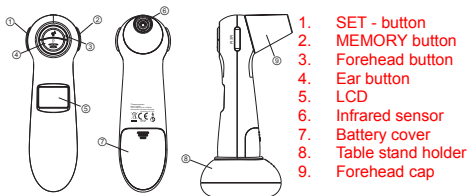
5 DISPLAY

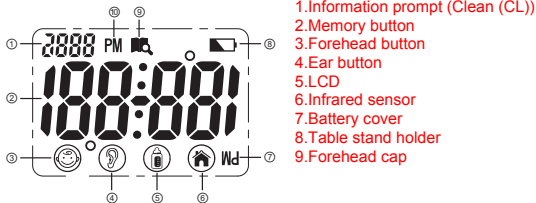6 BATTERY INSTALLATION

1. Remove the battery cover.
2. When the batteries are pre-installed. Pull out the insulator.
3. If it is not pre-installed, insert a CR2032 - 3V - Lithium battery with the positive pole (+) facing upwards.
4. Reposition the battery cover.

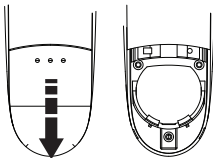IMPORTANT NOTE

- Replace the battery when '  ' is displayed.
- Take out the battery if you do not use it for a long time.
- You have to dispose of the batteries in an environment friendly manner according to your country regulations.

7 BODY TEMPERATURE

IMPORTANT NOTE

- It is recommended that you take three temperatures under the following circumstances.
 - Infants in the first 3 months of life.
 - Children under three years old who have a condition such as compromised immune system and for whom the presence or absence of fever is critical.
 - When you are first learning to use the thermometer until you are comfortable with the technique and are gaining consistent readings
- Clean the probe cover to ensure accurate measurement.
- Wait about 30 seconds before measuring the same person to avoid excessive cooling of the skin.
- As with any type of thermometer, slight temperature variations (+/- 0.2 - 0.3°C) can happen if consecutive measurements are taken.
- If the results are different, use the highest reading.
- **In measuring mode and memory mode, the LDCD is normal display; but in waiting mode and setting mode, the LCD is reverse display**
- **In waiting mode, the LCD circularly displays time, date and room temperature.**


7.1 HOW TO MEASURE FOREHEAD TEMPERATURE

1. Put on forehead cap before use. Place the thermometer in your hand as shown with your thumb on the [FOREHEAD] button.

- Do not press the button yet.
- Make sure the probe is clean before measuring

2. Hold the thermometer so that the forehead cap flush on the center of the forehead, midway between eyebrow and hairline. the forehead cap touch the forehead, press the [FOREHEAD] button to start the measurement

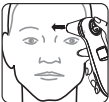3. Keeping the thermometer flat until you hear temperature broadcast or suggestive sound, and then the backlight will be light according to temperature. The measurement is complete and the thermometer may be removed from the head.

Read the temperature on the display.

7.2 HOW TO MEASURE EAR TEMPERATURE

1. Remove forehead cap before use. Place the thermometer in your hand as shown with your thumb or forefinger on the [EAR] button.

- Do not press the button yet.
- Make sure the probe is clean before measuring.

2. Insert the probe very gently and slowly into the ear canal. Hold the thermometer so that the probe faces straight in the direction of the eardrum. It is important to point the probe window of the thermometer toward the eardrum and its periphery inside the ear.

- Hold the outer ear and gently pull it toward the rear of the head to straighten the ear canal.

- Always insert the thermometer into the same ear with the same direction and depth.

3. Press the [EAR] button to start the measurement. Keeping the probe flat until you hear temperature broadcast or suggestive sound, and then the backlight will be light according to temperature.

4. Remove the thermometer from the ear canal. The LCD displays the measured temperature and the ear sign flash about 6 seconds. You can measure again after the ear sign stop flashing.

7.3 HOW TO MEASURE OBJECT TEMPERATURE

1. Remove forehead cap before use. Place the thermometer in your hand as shown with your fingers on the [EAR] and [FOREHEAD] button. Press the [EAR] and [FOREHEAD] button then the thermometer enter to object temperature measuring mode, the LCD display like left figure. Make sure don't use this mode to measure body temperature.

2. Hold the thermometer so that the probe position at a distance within 1cm from the liquid and /or surface whose temperature you want to measure.

Press the [EAR] button to start the measurement.

Remove the thermometer from the object. The LCD displays the measured temperature. You can press the [FOREHEAD] AND [EAR] button again at the same time or press [MEN] button return to waiting mode. The thermometer will automatically return to waiting mode after 1minute if you forgot to do so.

7.4 TEMPERATURE TAKING HINTS

1. A temperature measurement taken in the right ear may differ from a measurement taken in the left ear. Therefore, always take the temperature in the same ear.

2. External factors may influence ear and forehead temperatures, particularly when an individual has:

- been lying on one ear or the other
- had their ears covered
- been exposed to very hot or very cold temperatures
- been recently swimming or bathing.

In these cases, remove the individual from the situation and waiting for 30 minutes prior to taking a temperature.

3. For persons wearing hearing aids or ear plugs, remove the device and wait 30 minutes prior to taking a temperature.

4. Use the untreated ear if ear drops or other ear medications have been placed in the ear canal.

5. With object measurement, the current surface temperature of the object is displayed. It can be different from its internal temperature, especially if the surface is exposed to direct sunlight or a draught.

7.5 TEMPERATURE BACKLIGHT INDICATION

There are two temperature indication in this unit, the relevant backlight will light after human body temperature measurement finished. The detail range refer to following:

Blue backlight: below 37.9 °C (100.3.)

Red backlight: 38.0 °C ~42.9 °C (100.4.~109.3.)

7.6 BACKLIGHT FUNCATION

Press any button to turn on the backlight. The backlight will turn off automatically after 5 seconds or in temperature measuring mode.

7.7 CHECK THE MEMORY

There are a total of 9 memory data for recording ear and forehead temperature measurements.

The current measurement is always stored in the last storage space. When all storage spaces have been occupied, the oldest measurement is always deleted from the memory.

Press [MEM] button to check the memory number, temperature and date like right figures.

Press [MEM] button circulation to check temperature of the memory.

Press and hold [MEM] button to check the memory data fleetly.

Press [EAR] and [FOREHEAD] measuring button at the same time, the device back to waiting mode.

7.8 DATE AND TIME SETTING

Press and hold [SET] button about 2 seconds to enter into setting mode under the waiting mode.

Press [SET] button to switch the settings of hour mode (12hr or 24hr, in afternoon, "PM" will display on LCD when you select 12hr mode), hour, minute, year (2001~2099), month, day and sleep mode (ON or OFF).

When the sleep mode is active, the thermometer will turn off automatically about 1min or 3min(which is factory default); when the sleep mode is deactivated, the thermometer remains power on.

For save power, we recommended to set sleep mode in ON.



Press [MEM] button to increase the parameter which you selected.

Press [EAR] and [FOREHEAD] measuring button at the same time, the device back to waiting mode.

7.9 CELSIUS/FAHRENHEIT OPTION

Press and hold [SET] button to switch between Celsius(°C) and Fahrenheit(°F) when the thermometer in measuring mode or memory mode.

8 TROUBLESHOOTING

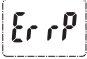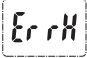Troubles	Checklist	Countermeasures
<p>No response / automatic reset when pull out insulator.</p>	<p>Battery used up</p> <p>Battery in wrong polarity</p> <p>Poor battery contact</p>	<p>Change new battery</p> <p>Take out battery, reinsert battery correctly.</p>
	<p>The measured temperature is lower than 32.0°C/ 89.6°F or higher than 42.9°C/ 109.3°F please check the operation method.</p>	<p>Follow steps in the user manual for proper measurement.</p>
	<p>Hardware problem</p>	<p>Contact the manufacturer or distributor</p>
		
		

The measurement is not accurate or if there is any doubt of the measure result

Please check if the infrared sensor is clean or not.

Clean the infrared sensor with cotton swab according to use the user manual.

Please check the measuring way is correct or not.

Ensure you have read the manual and know how to use the thermometer properly.

Please check if the thermometer and patient gets stabilized in the room for 30 minutes at least.

Please keep the thermometer and patient in the measuring room at least 30 minutes before using.

The measurement is not accurate or if there is any doubt of the measure result.

Are you using the thermometer indoor.

Please take the measurement indoor.

Please check if you hold the thermometer in your hand for too long this could affect the accuracy.

Put the thermometer on the table in the room where the measurement is taking place and let it cool down first.

Please check if the low battery icon is on the display

Replace the battery.

9 FEATURES

Model	TH-4654 + TH-4655
Temperature Sensor type	Infrared
Measuring temperature range	32 °C ~ 42,9 °C (89,6 °F ~ 109,3 °F)
Room temperature	5 °C ~ 59,9°C (41 °F ~ 139,8 °F)
Display accuracy	0,2°C or 0,4°F
Measurement accuracy	35,5°C ~ 42°C - (+/- 0.2°C) 32° C to 35.4 °C - (+/- 0.3 °C) 42,1C to 42.9 C° - (+/- 0.3 °C) 95,9°F ~ 107,6°F -(+/- 0,4°F) 89,6°F ~ 95,8°F - (+/- 0,5°F) 107,7°F ~ 109,3°F - (+/- 0,5°F)
Room temperature:	+/- 1°C
Food temperature:	+/-1°C
Operating ambient temperature range	15°C ~ 35°C (59°F ~ 95°F)
Operating ambient humidity range	30 ~ 80% RH
Storage environment temperature range	-10°C ~ 55°C (14,0°F ~ 131,0°F)
Storage environment humidity range	30 ~ 85% RH
Memory	Last 30 measurements
Battery Type	1 x CR2032 - DC 3V - Lithium
Automatic switch-off	3 minutes after the last button pressed in time mode
Weight	65g (with battery)
Classification	<ul style="list-style-type: none"> • Type B • Class IIa • Internally powered equipment • IPX0 • Not suitable for use in presence of flammable anesthetic mixture with air or with oxygen or nitrous oxide • Continuous operation with short-time loading

10 CE MARK

The CE symbol on the unit, user guide and giftbox indicates that the unit complies with the essential requirements of the medical directive 93/42/EEC. The declaration of conformity can be found on:

www.tristar.eu

11 CLEANING

- For reminding user to clean the probe tip, this thermometer has a special design—cleaning tip. When you are not use the thermometer has a long time or measure temperature a few times, the blue back light flashing, cleaning symbol(CL) display on LCD and DiDi warning sound at the same time to remind the user to keep clean the probe tip.
- The thermometer will return to original mode after press any button, and the cleaning indicator stop.
- To clean the tip of the probe, hold the thermometer downwards so that no liquid can enter. Use a cotton swab moistened with alcohol to wipe gently and carefully on the surface of the tip of the probe. Allow an hour drying time before operating the thermometer again.
- Clean a polluted infrared sensor with a soft wet cloth in a gentle manner. Cleaning with toilet paper or paper towel may scratch the infrared sensor, causing inaccuracy.

12 CARE AND MAINTENANCE

- Do not shock, vibrate or drop the thermometer.
- Do not allow any liquid in the thermometer.
- Do not disassemble, repair or remodel the thermometer.
- Keep out of reach of children.
- Clean the main body of the thermometer with a soft and dry cloth.
- Do not use benzene, gasoline, thinner or any other cleaner containing abrasive agents.

13 DISPOSAL OF THE DEVICE (ENVIRONMENT)

At the end of the product lifecycle, you should not throw this product into the normal household garbage but bring the product to a collection point for the recycling of electrical and electronic equipments. The symbol on the product, user guide and/or box indicate this.

Some of the product materials can be re-used if you bring them to a recycling point. By re-using some parts or raw materials from used products you make an important contribution to the protection of the environment. Please contact your local authorities in case you need more information on the collection points in your area.

Batteries must be removed before disposing of the device.

Dispose of the batteries in an environmental manner according to your country regulations.

14 TOPCOM WARRANTY

14.1 WARRANTY PERIOD

The Topcom units have a 24-month warranty period. The warranty period starts on the day the new unit is purchased. There is no warranty on standard or rechargeable batteries (AA/AAA type).

Consumables or defects causing a negligible effect on operation or value of the equipment are not covered.

The warranty has to be proven by presentation of the original or copy of the purchase receipt, on which the date of purchase and the unit-model are indicated.

14.2 WARRANTY EXCLUSIONS

Damage or defects caused by incorrect treatment or operation and damage resulting from use of non-original parts or accessories are not covered by the warranty.

The warranty does not cover damage caused by outside factors, such as lightning, water and fire, nor any damage caused during transportation. No warranty can be claimed if the serial number on the units has been changed, removed or rendered illegible.

Any warranty claims will be invalid if the unit has been repaired, altered or modified by the buyer.

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Shenzhen Donglixin Technology Co., LTD.
Shanghai International Holding Corp. GmbH (Europe)