MANUAL (p. 2)

ANLEITUNG (S. 6)

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MANUALE (p. 21)

KÄYTTÖOHJE (s. 26)

BRUKSANVISNING (s. 30)

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
INFRA-RED EAR THERMOMETER

How to use:

Measuring human body temperature in the ear canal

1. Press the on/recall button to turn on the thermometer, a beep sound is heard and the LCD displays as Fig A in one second. When the LCD displays as Fig B, now it's ready to measure.
2. Straighten the ear canal by pulling the outer ear up and back to give a clear view of the eardrum.
 - For children under 1 year, pull the ear straight back.
 - For children ages 1 year to adult, pull the ear up and back.

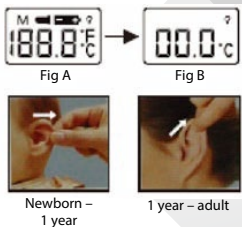
(Caution: Many pediatricians suggest use of ear thermometers for children older than 6 months).

3. While pulling the outer ear, insert the probe snugly into the ear canal (best to insert deep), press the "scan" button. Release it, when you hear a beep sound. A beep indicates measurement is complete. Remove the thermometer from the ear. The display shows the measure temperature.
4. Measuring again: If measuring again, please unclinch the "scan" button and wait for the "?" symbol to display, then press "scan" button to measure again.

Note: The time between each reading should be not less than 1 min.

5. Memory search: Press the "on/recall" button, the last ten memories (No. 9–No. 0) will take turns to display on the LCD as Fig C and Fig D.
6. The LCD will display "Er 1" as Fig E and cannot display temperature when the ambient temperature is of the range of 10 °C to 35 °C.

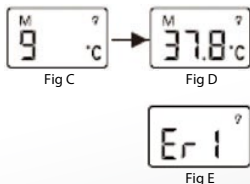
Note: Before measuring, the thermometer shall be stabilized at the operating ambient condition for a minimum of 30 min.

7. Shut off: The thermometer shall automatically shut off without

manipulating in one minute.

Note:

- In order to avoid spreading germs, please use a new disposable probe cover for each measurement.
- Check to ensure probe cover is fitted on firmly before use (please see the diagrams right). If the probe cover is broken or discarded, apply a new one immediately.

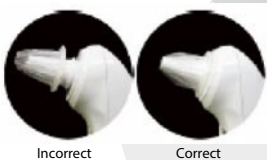Incorrect

Correct

Cleaning and storage:

1. Store thermometer and the disposable probe cover in a dry location free from dust and contamination and away from direct sunlight. The ambient temperature at the storage location should remain fairly constant and within the range of -10 °C to 55 °C.
2. Use an alcohol swab or cotton swab moistened with alcohol (70%) to clean the thermometer casing and the measuring probe. Ensure that no liquid enters the interior of the thermometer.
3. Never use abrasive cleaning agents, thinners or gasoline for cleaning and never immerse the instrument in water or other cleaning liquids. Take care not to scratch the surface of the probe membrane or display.

Replacing the batteries:

1. When voltage of the batteries are low, the LCD will display “” symbol, please replace two new batteries in the unit. The thermometer cannot work accurately under the condition of low voltage.
2. The thermometer is supplied with two 1.5 V AAA alkaline batteries. Insert new batteries when the low voltage symbol appears on the LCD.
3. Place two new batteries according to the + and -.
4. Please take out the batteries to avoid battery leaking if unit not used for over six months.

Specification:

The standard we adopted is EN12470-5:2003 clinical thermometers – Part 5: performance of infrared ear thermometers (with maximum device).

This device complies with the requirements of IEC 60601-1:1998, and IEC 60601-1-2:2001

Range of displayed temperature: 34 °C–44 °C

Display L°C when the temperature is under 34 °C

Display H°C when the temperature is over 44 °C

Operating ambient:


Temperature: 16 °C–35 °C

Humidity: < 80% RH

Storage and transport condition:

Temperature: -10 °C–55 °C

Humidity: 30–90% RH

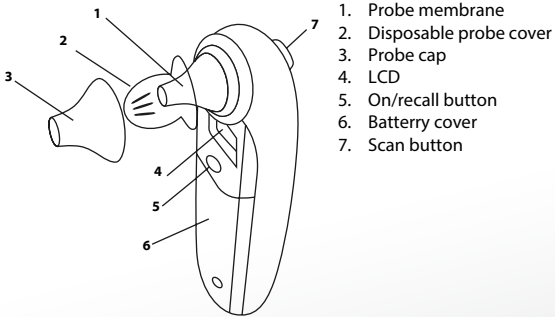Display resolution temperature range: 0.1 °C**Accuracy: +/-0.2 °C (35.5 °C–42.0 °C)****Display:** Liquid crystal display, 4 digits**Power consumption:** 0.6 million watt in measurement mode**Battery:** Two 1.5 V alkaline batteries (AAA)

 **Cautions**

1. Before measuring, the thermometer shall be stabilized at least 30 min under operating condition.
2. New disposable probe cover required for each measurement, otherwise inaccurate readings will occur.
3. Please keep your ear canal clear, or the measuring result shall be inaccurate.
4. The symbol “?” on the LCD means you can measure in the ear.
5. Please don't scratch the probe membrane, or the thermometer shall lose the efficacy.
6. If the LCD has no changes under pressing any button, please take the batteries out and put in again.
7. Please don't use the thermometer if your ear canal is inflamed.
8. “⚡”: Degree of protection against electric shock is type BF applied part.
9. “♻️”: The disposal of battery and device shall comply with the local environment requirements.

**Maintenance:**

Clean only with a dry cloth. Do not use cleaning solvents or abrasives.

Warranty:

No guarantee or liability can be accepted for any changes and modifications of the product or damage caused due to incorrect use of this product.

General:

Designs and specifications are subject to change without notice. All logos, brands and product names are trademarks or registered trademarks of their respective holders and are hereby recognized as such.

Attention:

This product is marked with this symbol. It means that used electrical and electronic products should not be mixed with general household waste. There is a separate collections system for these products.

This product has been manufactured and supplied in compliance with all relevant regulations and directives, valid for all member states of the European Union. It also complies to all applicable specifications and regulations in the country of sales.

Formal documentation is available upon request. This includes, but is not limited to: Declaration of Conformity (and product identity), Material Safety Data Sheet, product test report.

Please contact our customer service desk for support:

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NEDIS B.V., De Tweeling 28, 5215 MC 's-Hertogenbosch, NEDERLANDENE