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1 safety instructions

Important safety instructions - read carefully and keep for later use!

WARNING

- Failure to observe the information provided below can lead to personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. Therefore, observe this safety information and pass these instructions together with the product.
- This heating pad must not be used by people who are not sensitive to heat or by vulnerable people who may not be able to react to overheating.
- This heating pad may not be used by very young children (0-3 years old) due to their inability to react to overheating.
- This heating pad may not be used by young children (3-8 years old) unless the controller has been pre-set by a parent or guardian and unless the child has been adequately instructed how to operate the heating pad safely.
- This heating pad may be used by children over the age of 8 and by people with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the heating pad.
- Cleaning and user maintenance shall not be made by children without supervision.
- This heating pad is not intended for use in hospitals.
- Do not insert pins.
- Do not use folded or bunched up.
- Do not use if wet.
- This heating pad must be used with the type of controller that is marked on the rating label of the heating pad, only.
- The electrical and magnetic fields emitted by this heating pad may interfere with the function of a pacemaker. However, they are still well within the limits: electrical field strength: max. 5000 V/m, magnetic

field strength: max. 80 A/m, magnetic flux density: max. 0.1 millitesla.
Please consult your doctor and the manufacturer of your pacemaker before using this heating pad.

- Do not pull, twist or make sharp bends in the cables.
- This heating pad must be examined frequently for signs of wear or damage. If there are such signs, if the heating pad has been mis-used or does not work, return it to the supplier before switching it on again.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similar qualified persons in order to avoid a hazard.
- While this heating pad is switched on
 - never place any objects on it (e.g. suitcase or washing basket)
 - never place any heat source on it (e.g. hot water bottle, heating pad or similar).
- When in use, the electronic components in the control of the heating pad make the control slightly warm. The control should therefore never be covered or placed on the heating pad when it is in operation.
- If you should have any questions about using our devices, please contact Customer Services.
- Always observe the notes on operation (section 4), cleaning and maintenance (section 5) and storage (section 6).

2 appliance description

HNK115

Elektric heating pad

Dimensions: 47 x 34 cm

6 illuminated settings

Turbo heating


Automatic switch-off after 90 minutes

Detachable cord

Machine washable, 30°C

100 Watt

Material: fleece

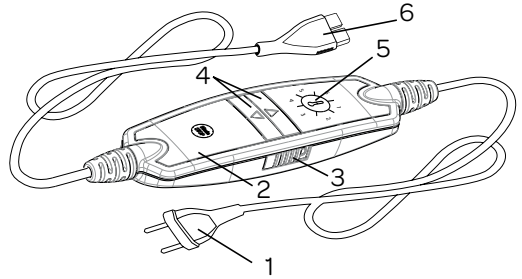1. Supply plug
2. Controller
3. Sliding switch (ON = I/ OFF = 0)
4. Buttons for setting temperature
5. Illuminated display for temperature settings
6. Plug-in coupling

3 intended use

CAUTION

This heating pad is only designed to warm up the human body.

4 operation

4.1 Safety

CAUTION

This heating pad is fitted with the INVENTUM SAFETY SYSTEM (ISS). This electronic sensor technology prevents the entire surface of the heating pad from overheating by automatically switching off in the event of a fault. If the safety system switches off the heating pad automatically as a result of a fault, the display for the temperature settings is no longer illuminated even when the heating pad is switched on.

Please note that, after a fault has occurred, the heating pad can no longer be used for safety reasons and must be sent to the service address indicated.


4.2 Initial use

CAUTION

Ensure that the heating pad will not bunch up or become folded during use.

- First, connect the controller to the heating pad by inserting the plug-in coupling.
- Then plug the mains plug into the socket.

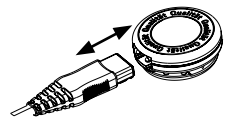4.3 Switching on

Push the sliding switch on the right side of the controller to the setting I (ON). When the switch is on, the temperature settings display is illuminated.

4.4 Setting the temperature

To increase the temperature, press the Δ button. To reduce the temperature, press the ∇ button.

Setting 1: minimum heat
Setting 2-5: individual heat setting
Setting 6: maximum heat

NOTE

The fastest way to warm up the heating pad is to initially set the highest temperature setting of the controller, first.

NOTE

This heating pad has a turbo function, which allows the heating pad to warm up faster in the first 10 minutes.

CAUTION

If the heating pad is being used over several hours it is recommended to use the lowest setting (1) of the control to avoid overheating the heated body part, which may lead to burns to the skin.

4.5 Automatic switch-off

This heating pad is equipped with an automatic switch-off function. This turns off the heat supply approx. 90 minutes after the initial use of the heating pad. A part of the displayed temperature settings on the controller then begins to flash. So that the heating pad can be switched back on, the side sliding switch (3) must first be set to setting "0" (OFF). After about 5 seconds it is possible to switch it on again.

4.6 Switching off

To switch the heating pad off, set the sliding switch for ON/OFF to setting 0 (OFF) position. The temperature settings are then no longer illuminated.

NOTE

If the heating pad is not in use, switch the side sliding switch for ON/OFF to the OFF (0) position and remove the mains plug from the socket. Then disconnect the controller from the heating pad by unplugging the plug-in coupling.

5 cleaning & maintenance

WARNING

Before cleaning, always remove the mains plug from the socket first. Then disconnect the controller from the heating pad by unplugging the plug-in coupling. Otherwise there is a risk of electric shock.

CAUTION

• The controller must never come into contact with water or other liquids. Otherwise, the controller may be damaged.

To clean the controller, use a dry, lint-free cloth. Do not use any chemical or abrasive cleaning agents.

The textile cover can be cleaned in accordance with the symbols on the label and must be removed from the heating pad prior to cleaning. Small marks on the heating pad can be removed with a damp cloth or sponge, or if necessary, with a little liquid detergent for delicate laundry.

CAUTION

• Please note the heating pad must not be dry-cleaned, wrung out, machine dried, pressed or ironed.

In case of heavier soiling, this heating pad can be washed in a washing machine. Set the washing machine to a particularly gentle, 30°C cycle (woolens programme). Use a wool laundry detergent and follow the manufacturer's instructions.

CAUTION

- Please note that excessive washing will wear out the heating pad. For this reason, the heating pad should only be machine-washed a maximum of 10 times throughout its entire life.

Immediately after washing, reshape the heating pad to its original dimensions while it is still damp and spread it out flat over a clothes drying rack to dry.

CAUTION

- Do not attach the heating pad to washing lines with washing pegs or similar. Otherwise the heating pad may be damaged.
- Do not reconnect the switch to the heating pad until the plug-in connection and the heating pad are completely dry. Otherwise the heating pad may be damaged.

WARNING

- Never switch the heating pad on to dry it! Otherwise there is a risk of electric shock.

6 storage

If you do not plan to use the heating pad for a long period of time, we recommend that you store it in the original packaging.

CAUTION

Please allow the heating pad to cool down first. Otherwise, the heating pad may be damaged. To avoid sharp folds in the heating pad, do not place any objects on top of it while it is being stored.

7 disposal

Legislation requires that all electrical and electronic equipment must be collected for reuse and recycling. Electrical and electronic equipment marked with the symbol indicating separate collection of such equipment must be returned to a municipal waste collection point.

8 problems & solutions

Problem	Cause	Remedy
The temperature settings are not illuminated while - the controller is fully connected to the heating pad; - the supply plug is connected to the mains outlet with correct function; - the sliding switch on the controller is set to setting I (ON).	The safety system switched off the heating pad irreversibly.	Send the heating pad together with the controller to the service department.

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guarantee

- This product is guaranteed for a period of 2 years from the date of purchase against mechanical and electrical defect.
- Guarantee can only be given if proof of purchase (e.g. till receipt) can be shown. Therefore keep it in a safe place. The guarantee period starts at date of purchase.
- Guarantee is not valid in case of:
 - normal wear;
 - damage through misuse use or when the instructions were not followed correctly;
 - lack of maintenance;
 - not following the operation- and safety instructions;
 - incompetent montage or repairs made by yourself or others;
 - using accessories and spare parts which are not suitable for the appliance;for assessment by our technical department.

- The guarantee is not valid for:
 - normal wear;
 - external cables;
 - (interchangeable) signal- and controllampsfor assessment by our technical department.

after sales service

Return to your supplier or send the appliance, including a copy proof of purchase and description of the complaint, stamped to:

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When the guarantee is no longer valid or when the period of guarantee has passed, we will charge any costs that flow out of an assessment and/or repairing.

On request we will supply a statement before repairing. Without requesting a statement, we assume that you, by sending the appliance, automatically order to repair the appliance.

A guarantee period of two months is valid for repairs done by our technical department and the used spare parts, counting from the repair date. When a malfunction occurs within this period which directly flows out of the repairs and the used spare parts, the consumer has the right of repairs without any costs.

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Inventum Huishoudelijke Apparaten B.V.
Apparaten B.V.**