Silk'n®

Beautiful Technology

OPTIFLOSS





ND-DC01

EN
DE
FR
NL
ES
PT
GR
TR
IT
NO
SV
DA

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DE USUARIO
- PT MANUAL DO UTILIZADOR

- GR ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ
- TR KULLANIM KILAVUZU
- IT MANUALE D'USO
- NO BRUKERHÅNDBOK
- SV ANVÄNDARMANUAL
- DA BRUGERMANUAL

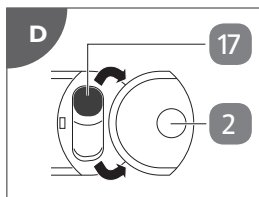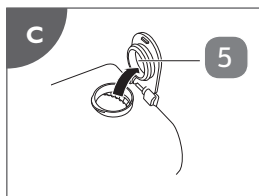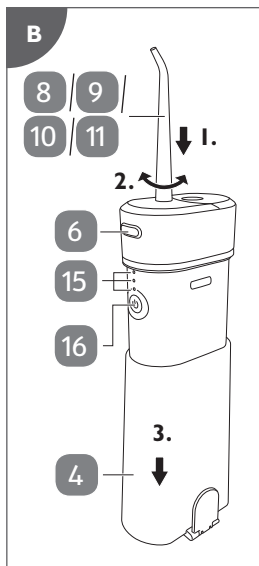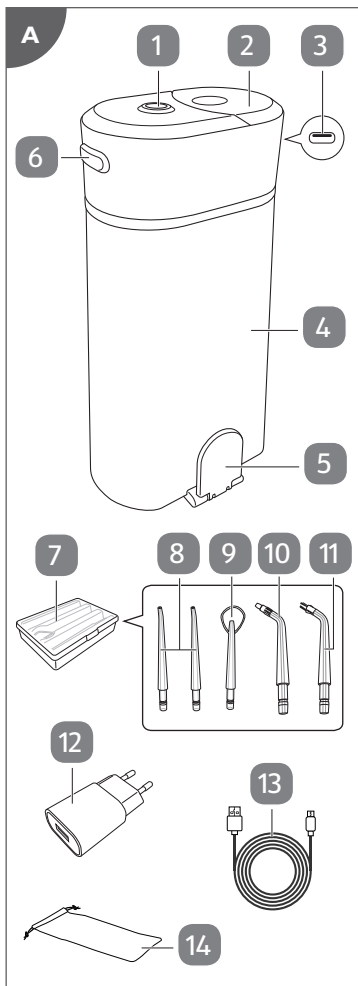Legal notice

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

SILK'N and the SILK'N logo are registered trademarks of Home Skinovations Ltd.

Invention Works B.V.

Donk IB, 2991 LE Barendrecht, The Netherlands

www.silkn.eu – info@silkn.eu

Table of contents

1	Package contents/device parts	5
2	General information	5
2.1	Reading and storing the user manual	5
2.2	Explanation of symbols	5
3	Safety	7
3.1	Proper use	7
3.2	Safety instructions	7
4	Checking the product and package contents	11
5	Charging	11
6	Operation	12
6.1	Oral cleaning	12
6.2	After every use	13
6.3	Removing or replacing the tip	14
7	Cleaning	14
7.1	Main unit	14
7.2	Water tank	14
7.3	Tips	15
8	Storage	15
9	Troubleshooting	15
10	Technical data	16
10.1	General	16
10.2	USB wall adapter	16
11	Disposal	17
11.1	Disposing of the packaging	17
11.2	Disposing of the product	18
12	Warranty	18
13	Customer service	19

I Package contents/device parts

- | | |
|-----------------------------|----------------------------|
| 1 Connection hole | 9 Tongue cleaner tip |
| 2 Storage compartment cover | 10 Orthodontic tip |
| 3 Charging port | 11 Plaque remover tip |
| 4 Water tank | 12 USB wall adapter |
| 5 Water tank cover | 13 USB cable |
| 6 Tip release button | 14 Travel pouch |
| 7 Storage case for tips | 15 Indicator light, 3× |
| 8 Standard tip, 2× | 16 Operation button (⏻) |
| | 17 Tip storage compartment |

2 General information

2.1 Reading and storing the user manual

This user manual accompanies this OptiFloss portable dental waterflosser (hereafter referred to as the "product"), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product. The user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

This symbol identifies voltage and the current rating for direct current.

The word "GS" stands for tested safety. Products marked with this symbol meet the requirements of the German Product Safety Act (ProdSG).

IP20

The USB wall adapter is protected against solid objects up to 12 mm.

The international efficiency mark indicates that the power supply meets the level VI requirements (USB wall adapter).

The USB wall adapter is equipped with a short-circuit proof safety transformer.

The USB wall adapter has a switching power supply.

This symbol indicates the polarity of the pins in the USB wall adapter.

This symbol identifies electrical devices that belong to the protection class II.

ta

Ambient temperature (USB wall adapter)

3 Safety

3.1 Proper use

The product is exclusively designed to remove plaque and food residues, and to maintain your oral hygiene. It is only intended for private use and not suitable for commercial purposes.

Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer cannot be held liable for damages or injury incurred through improper or incorrect use.

3.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- Only connect the product if the line voltage of the socket corresponds to the data on the rating plate of the USB wall adapter.
- Only connect the USB wall adapter to an easily accessible socket so that you can quickly disconnect it from the mains supply in the event of a problem.
- The supplied USB wall adapter can only be connected to the product. Do not use it for other purposes.
- Do not use the product, USB cable or USB wall adapter if they are damaged or defective. The manufacturer or customer service must replace the damaged or defective parts to avoid risks.
- Do not open the housing; instead, have a qualified professional perform repairs. Contact the customer service team for this. Liability and warranty claims are waived

in the event of repairs performed by the user, improper connection or incorrect operation.

- Only parts that comply with the original device data may be used for repairs. This product contains electrical and mechanical parts which are essential for providing protection against sources of danger.
- Do not immerse the USB cable or USB wall adapter in water or other liquids.
- Do not leave the product submerged in water (max. 1 m) or under running water for more than 30 minutes.
- Always switch the product off when you are not using it, before changing the tip, when you intend to clean it, or in the event of a malfunction. Always unplug the USB wall adapter when you are not charging the product.
- Never touch the USB cable and USB wall adapter with wet or damp hands.
- Never insert any objects into the housing.
- Do not pull the USB wall adapter out of the socket by the USB cable; instead, always pull it out by the adapter itself. Never move, pull, or carry the product by the USB cable.
- Keep the product, USB cable and USB wall adapter away from open flames and hot surfaces.

- Lay the USB cable so that it does not pose a tripping hazard.
- Do not kink the USB cable or wrap it around the product. Do not lay the USB cable over sharp edges.
- Only use the product indoors.

WARNING! RISK OF INJURY!

- The product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the product in a safe way and understand the hazards involved. Children must not play with the product. Cleaning and user maintenance must not be performed by children unless they are older than 8 years and supervised.
- Keep children under the age of 8 away from the product, USB cable and USB wall adapter.
- Make sure that children do not play with the plastic wrapping. They may get caught in it when playing and suffocate.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

- Use only cold or warm water (not exceeding 40 °C) with the product.
- Do not leave the product unattended while in use.
- Do not direct the water pulses under the tongue or into the ear, nose, eye, or other sensitive areas.
- Do not press the tip too firmly against your teeth or gums.
- Consult your dentist before you use this product if you have had recent oral or gum surgery, suffer from an ongoing dental condition, had extensive dental work, have suspected issues with your dental work, or are not sure whether the product is suitable for you (e.g. using a pacemaker).
- Stop using the product and contact a physician/dentist if discomfort or pain is experienced.

NOTICE! RISK OF DAMAGE!

- Keep the product away from open fire, hot surfaces, high temperatures or the effects of the weather (e.g. rain).
- Do not place the product on the edge of a surface to avoid dropping.
- Unplug this product before filling the water tank. Fill the water tank with water or mouthwash only. Do not overfill.

- Never fill the housing with liquid.
- Stop using the product if its plastic parts exhibit cracks or breaks or are deformed. Only replace damaged parts with corresponding original spare parts.

4 Checking the product and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
1. Take the product out of the packaging.
 2. Check to make sure that all parts are included (see **Fig. A**).
 3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

5 Charging

WARNING! RISK OF INJURY!

- Charge only with the supplied USB cable and USB wall adapter.
- Only charge the product at temperatures between 0 and 40 °C.

- Fully charge the product before first use and when the battery power is low.
- The top indicator light  flashes during operation when the battery power is low. Recharge the product.
- The product cannot be used during charging.

1. Plug the smaller end of the USB cable **13** into the charging port **3** of the product.
2. Plug the larger end of the USB cable into the USB port of the USB wall adapter **12**.

You may also connect the USB cable to a compatible USB power source (e.g. PC, laptop or powerbank) for charging.

3. Plug the USB wall adapter into a mains socket for charging.

All indicator lights flash when charging starts.

The bottom indicator light stops flashing and stays lit when the charging level reaches 60 %.

The middle indicator light stops flashing and stays lit when the charging level reaches 80 %.

All three indicator lights stay lit when the charging level reaches 98 %.

4. When charging is done, unplug the USB wall adapter from the mains socket and disconnect the USB cable from the product.

6 Operation

6.1 Oral cleaning

1. Insert a tip into the connection hole **1**. You will hear a click sound when the tip is fully inserted (see step 1 in **Fig. B**).

Use a standard tip **8** for general teeth and gum cleaning.

There are two standard tips (0.6 / 0.8 mm). The one with a smaller diameter offers higher water pressure.

Use the tongue cleaner tip **9** for tongue scraping and cleaning.

Use the orthodontic tip **10** if you have braces, or other hard-to-reach places.

Use the plaque remover tip **11** to remove plaque around the gums.

2. Rotate the tip so that it is facing a direction convenient for cleaning (see step 2 in **Fig. B**). The tip can rotate 360°.
3. Gently pull the water tank **4** downwards until it clicks into place (see step 3 in **Fig. B**).
4. Hold the product horizontally and flip open the water tank cover **5** (see **Fig. C**).
5. Fill the water tank with water (not exceeding 40 °C) and close the cover. Do not overfill the water tank.

6. Hold the product in an upright position with the top part of the standard tip pointing towards the cleaning area (teeth or gum). Keep your mouth slightly opened so that water can flow out.
7. Press the operation button (⏻) 16 once.
The product is now switched on and in the Normal mode for normal cleaning. The top indicator light 15 illuminates.
Water pulses come out from the tip once the product is switched on.
8. Press the operation button (⏻) again.
The product is now in the Soft mode for gentle cleaning. The middle indicator light illuminates.
9. Press the operation button (⏻) a third time.
The product is now in the Pulse mode for cleaning and massaging. The bottom indicator light illuminates.
In this mode, the product uses alternating strong and weak water pulses for a massage-like soothing effect.
10. Press the operation button (⏻) a fourth time.
The product is now in the Manual mode for manual cleaning. All three indicator lights illuminate. In this mode, the ejection of water pulses can be manually adjusted. Press the operation button (⏻) again and hold the button until it has reached the desired amount of water pulses.
11. With the operation mode selected, direct the water pulses slowly along the teeth and gum line at a right angle.
For tongue cleaning, use the cup on the tongue cleaner tip to gently scrap against the tongue surface from the back to the front.
12. The product automatically switches off after 90 seconds.
You can also press the operation button (⏻) a fifth time, or press and hold the operation button (⏻) to switch the product off any time.
All indicator lights turn off when the product is switched off.

6.2 After every use

1. Open the water tank cover 5, pour out the remaining water from the water tank 4 and then close the water tank cover.
2. Direct the tip to a sink and press the operation button (⏻) 16 to switch the product on for a few seconds. This will drain the remaining water inside the water tank and the internal parts of the product.

3. Wipe the product dry with a dry cloth.
4. Gently push the water tank upwards until it is fully closed.

6.3 Removing or replacing the tip

1. Make sure the product is switched off.
2. Press and hold the tip release button **6** and lift the tip **8/9/10/11** upwards to remove it from the connection hole **1**.
3. Release the tip release button.

For tip replacement, insert another tip into the connection hole. You will hear a click sound when the tip is fully inserted.

7 Cleaning

WARNING! RISK OF ELECTRIC SHOCK!

- Always unplug the product from mains power before cleaning.

NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners (e.g. vinegar or limescale removers), brushes with metal or nylon bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers and the like.
- Do not use hot water for cleaning.

7.1 Main unit

Use a damp cloth (if necessary, with liquid soap) to wipe off the main unit and then wipe dry with a dry cloth.

7.2 Water tank

1. Clean the inside of the water tank **4** with warm water. Remove any dirt at the water tank cover **5** with a damp cloth.
2. Pour out the water from the water tank. Always keep the water tank dry when not in use and before being stored.

7.3 Tips

1. Clean the tips **8**/**9**/**10**/**11** with warm water.
2. Wipe the tips dry with a dry cloth.

8 Storage

All parts must be completely dry before being stored. Make sure there is no water in the water tank **4**.

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product in the travel pouch **14** when travelling to avoid damage.
- Store the tips **8**/**9**/**10**/**11** in the storage case for tips **7**. One of the standard tips can also be stored in the tip storage compartment **17** when travelling to avoid damage. Push the storage compartment cover **2** sideways to open it (see **Fig. D**).

9 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not operate.	<ul style="list-style-type: none"> • Make sure the product is properly charged.
The product does not charge.	<ul style="list-style-type: none"> • Make sure the USB cable is securely connected between the product and USB wall adapter, and the USB wall adapter is properly connected to a mains socket.
Water pressure is very low.	<ul style="list-style-type: none"> • Make sure the product is properly charged. • Make sure the tip is not dirty or deformed.
Water does not come out from the tip.	<ul style="list-style-type: none"> • Make sure the water tank is not empty. • Make sure the product is not tilted. Keep the product in an upright position.

10 Technical data

10.1 General

Model:	ND-DC01
Rated voltage:	DC 5 V
Power:	Max. 8 W
Rechargeable battery:	1,500 mAh lithium, 3.7 V
Charging time:	3 hours
Operating time:	Up to 60 minutes / full charge
Automatic power off:	90 seconds
Water tank capacity:	180 ml
Noise level:	≤ 70 dB
Pump frequency:	2,100 pulses/minute (normal mode) 1,260 pulses/minute (soft mode) 2,100 or 1,260 pulses/minute (pulse mode)
Water pressure:	6.2 bar (normal mode) 3.8 bar (soft mode)
IP code:	IPX7 (protected against water immersion for up to 30 minutes at a depth of 1 m)
Weight:	230 g
Dimensions:	4.9 × 7.45 × 14.35 cm

10.2 USB wall adapter

Model number:	PS06H050K0500EU
Input:	100–240 V~ 50/60 Hz 0.25 A
Output:	5.0 V === 0.5 A, 2.5 W
Protection class:	II
IP code:	IP20 (protected against solid objects up to 12 mm)

Published information	Value and precision	Unit
Name or trademark of the manufacturer, commercial register number and address:	Shenzhen Flypower Technology Co., Ltd. Building A2, Haosan Linpokeng 2nd Industry Zone, Nanpu Road, Shajing, Bao'an, Shenzhen, Guangdong, China	
Model identifier:	PS06H050K0500EU	
Input voltage:	100–240	V
Input alternating current frequency:	50/60	Hz
Output voltage:	5.0	V
Output current:	0.5	A
Output power:	2.5	W
Average efficiency during operation:	73.61 (min.)	%
Efficiency at low load (10 %):	63.55 (min.)	%
Power consumption at zero load:	0.0396 (max.)	W

II Disposal

II.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

11.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in the household waste!

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or**

county. This ensures that old appliances are recycled in a professional manner and also rules out negative

consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries may not be disposed of with household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an

environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

Note: To remove the rechargeable battery, wrap a towel around the product and use a hammer to smash the housing. Make sure your eyes and hands are protected to avoid injury. Please note, the removal of the rechargeable battery is not reversible.

12 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

13 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: www.silkn.eu.

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Copyright 2021 © Home Skinovations Ltd. Alle rechten voorbehouden.
Copyright 2021 © Home Skinovations Ltd. Todos os direitos reservados.
Copyright 2021 © Home Skinovations Ltd. Alle rettigheder forbeholdes.
© 2021 Home Skinovations Ltd. All rights reserved. SILK'N and the SILK'N logo are registered trademarks of Home Skinovations Ltd.